

INTEGRAPH™

Multi-Function Chemical Integrator for steam sterilization at 132°C (270°F)

Catalog No. CI 101

(EN ISO 11140-1, Class 5)

- **VIVID PURPLE TO GREEN COLOR CHANGE**

- **MULTI-FUNCTION WITH THREE INK CHEMISTRIES**

Check Mark—process indicator
Long Bar—integrator
Short Bar—extended cycle

- **SPECIAL FEATURES**

Turns dark brown under superheated steam conditions.


Will not complete a color change if excess air is present.

Plastic lamination to avoid direct contact with sterilized items.


- **MEETS HIGH STANDARDS**

Released for market by the FDA, meets ISO standards as a true integrator.


UNEXPOSED

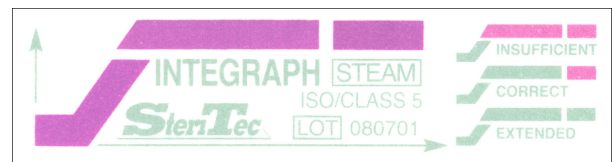

CORRECT COLOR CHANGE (NORMAL CYCLE)

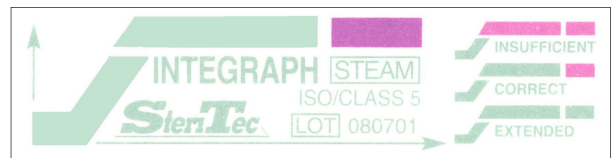

CORRECT COLOR CHANGE (EXTENDED CYCLE)

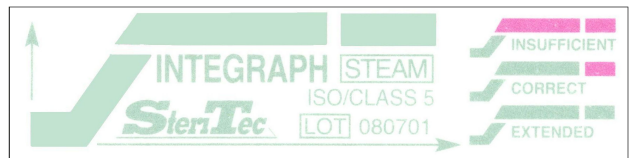

ORDERING INFORMATION

<u>PART #</u>	<u>PACKAGING</u>
CI 101	20BAG of 250/CASE

Steri-Tec

PRODUCTS, INC.

599 TOPEKA WAY • SUITE 400
CASTLE ROCK, CO 80109 • USA
(303) 660-4201 • Fax: (303) 660-4213
Website: <http://www.steritecproducts.com>
E-Mail: info@steritecproducts.com