

STERI-VAPEUR 18

18 Minute Integrator

Catalog No. CI 117

EN ISO 11140-1, Class 6

- **TWO COLOR CHANGE CHEMISTRIES**

The 4 MIN spot relates to normal cycle applications such as cloth, rubber goods, and instruments.

The 18 MIN spot relates to extended exposures.

- **SPECIAL FEATURES**

Plastic lamination protects sterilized items from potential ink transfer.

Requires both steam and heat, for color change to occur.


Stable chemistry. Colors will not reverse. Can be used as a permanent record.

Also indicates superheated and wet steam conditions.


UNEXPOSED



NORMAL CYCLE (4 MINUTES)

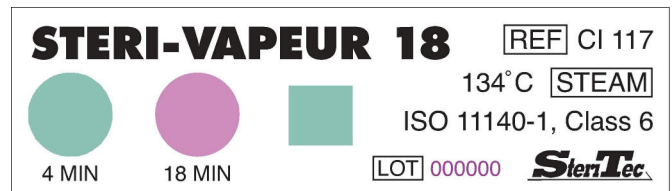

EXTENDED CYCLE (18 MINUTES)

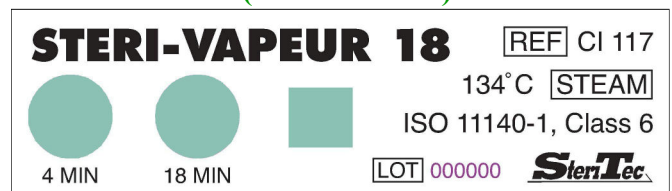

ORDERING INFORMATION

PART #	PACKAGING
CI 117	16 BAGS OF 200/CASE

Steri-Tec

PRODUCTS, INC.

599 TOPEKA WAY • SUITE 400
CASTLE ROCK, CO 80109 • USA
(303) 660-4201 • Fax: (303) 660-4213
Website: <http://www.steritecproducts.com>