

# EMU-GRAPH EMULATOR

Catalog No. CI 131

ANSI/AAMI/ISO 11140-1, Class 6

- **TWO COLOR CHANGE CHEMISTRIES**

The 2 MIN. bar relates to biological kill time at 134°C.

The 4 MIN. bar relates to sterility assurance level (SAL) of  $10^6$  which is twice biological kill as required in hospitals.

- **EMULATOR STATUS**

Completes color change within 6% of stated value as defined by ANSI/AAMI/ISO 11140-1.

- **SPECIAL FEATURES**

Plastic lamination protects sterilized items from potential ink transfer.

Requires both steam and heat, for color change to occur.

Stable chemistry. Colors will not reverse. Can be used as permanent record.

Also indicates superheated and wet steam conditions.


## ORDERING INFORMATION

PART #	PACKAGING
CI 131	20 BAGS OF 250/CASE

## BEFORE



## RELATES TO BIOLOGICAL KILL

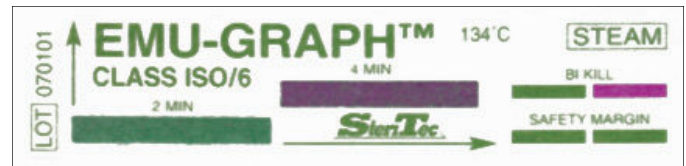

## SAFETY MARGIN



# SteriTec

PRODUCTS, INC.

599 TOPEKA WAY • SUITE 400

CASTLE ROCK, CO 80109 • USA

(303) 660-4201 • Fax: (303) 660-4213

Website: <http://www.steritecproducts.com>